
 <p>Advanced Accelerator Applications A Novartis Company</p>	<p><b>MAINTENANCE TRACKING TOOL</b></p> <p><b>PETTRACE800</b></p>	<p>Date:2024-06-03</p>
---	---	------------------------

<p><b>Country:</b> Italy</p>	<p><b>Site:</b> VEN</p>
<p><b>Intervention:</b></p>	<p>Programmed maintenance: UBM/CBM <input checked="" type="checkbox"/></p>
<p><b>Subsystems:</b></p>	<p></p>

**Vacuum leak test**

Seconds since push standby	Chamber pressure	Max. Chamber pressure
0	1.10E-7	1.80E-07
10	6.00E-7	1.00E-06
20	1.00E-6	1.50E-06
30	1.20E-6	1.90E-06
40	1.50E-6	2.30E-06
60	1.60E-6	3.00E-06
70	1.80E-6	3.30E-06
80	2.00E-6	3.60E-06
90	2.10E-6	3.90E-06
100	2.30E-6	4.20E-06
110	2.50E-6	4.60E-06
120	2.70E-6	4.90E-06

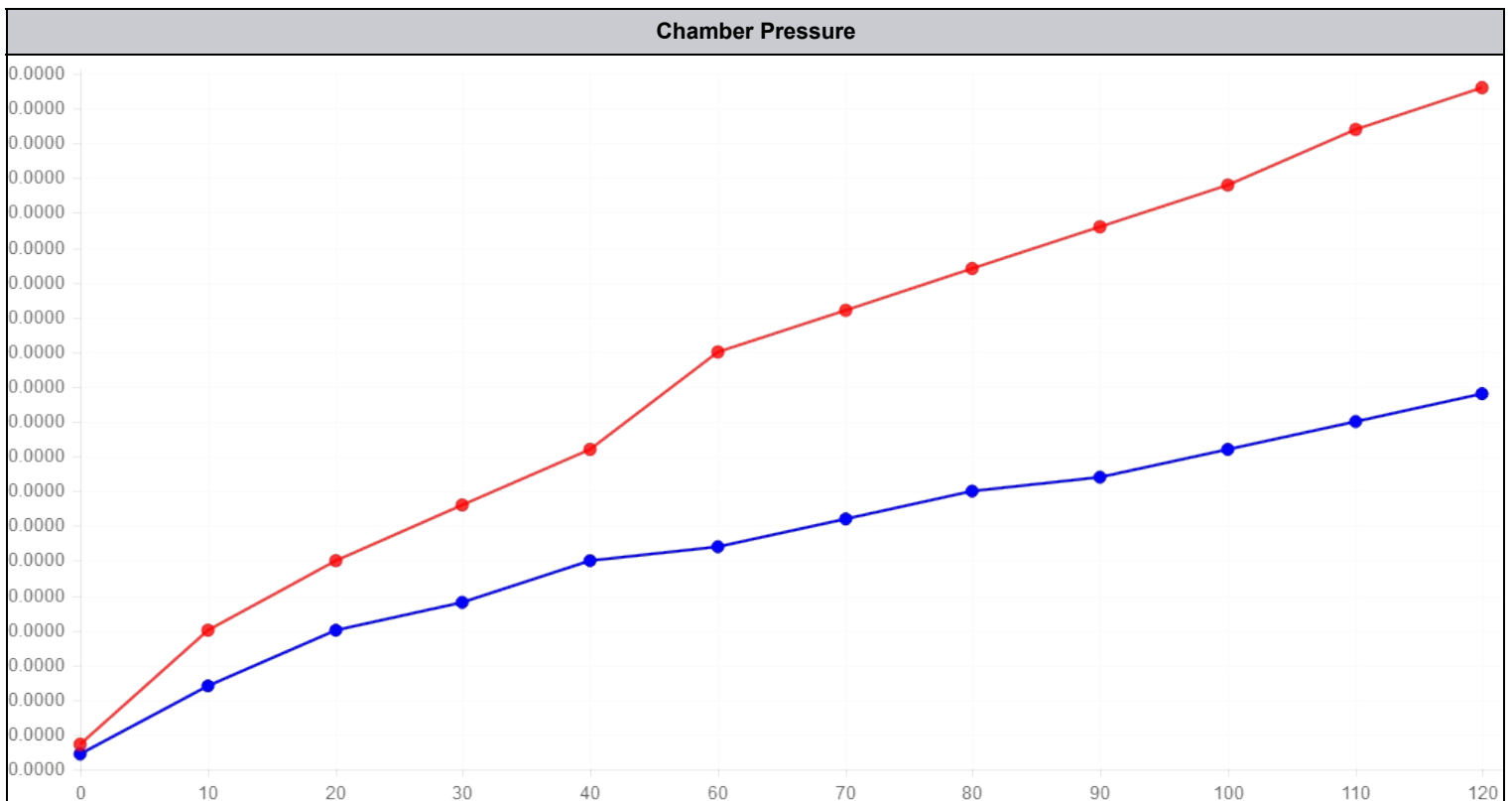

PETtrace800 O-Rings analysis

Pins

