



MAINTENANCE TRACKING TOOL

PETTRACE800

Date:2023-11-13


Country: Italy	Site: VEN
Intervention:	Single intervention <input checked="" type="checkbox"/>
Subsystems:	Chamber <input checked="" type="checkbox"/>

CHAMBER

Chamber Opening

Measure yoke play, adjust if needed: 2.0


Dose rate mapping (positions 1-9, [$\mu\text{Sv/h}$])

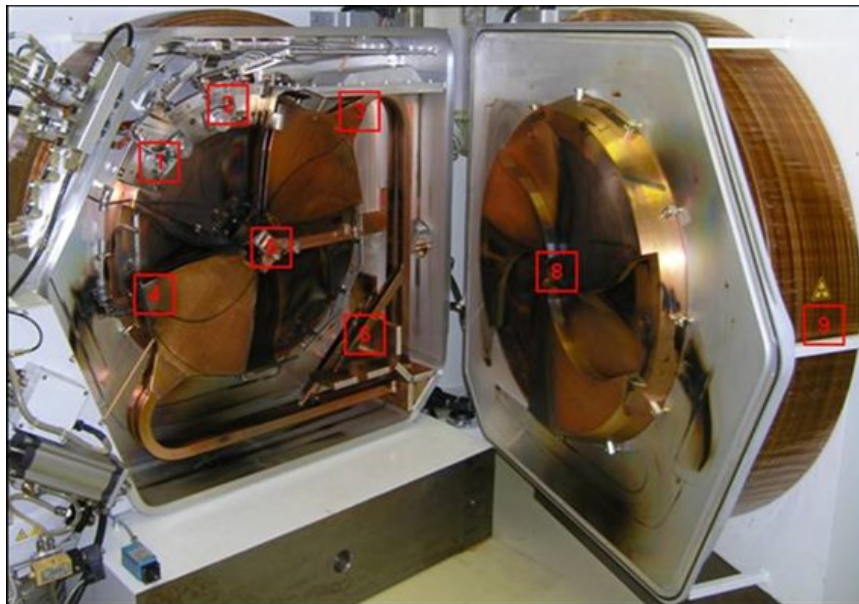

Position 1: At 36 cm from Extraction trolley	170
Position 2: At 36 cm from Carousel	170
Position 3: At 36 cm from Dee 2-stem junction	160
Position 4: At 36 cm from Deel upper corner	190
Position 5: At 36 cm from Central region	160
Position 6: At 36 cm from Stems coupler	140
Position 7: At contact with central region	190
Position 8: At 36 cm from magnet pole	110
Position 9: At contact of magnet coil	230

Photo documentation & visual inspection

Image



Beam exit valve tests

	<input checked="" type="checkbox"/>
--	-------------------------------------

